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## Piston Compressor Nebulizer "TeddyNeb"

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## Introduction

---

Congratulations on your purchase of the San Up Piston Compressor Nebulizer "Teddy-Neb". Thanks for trusting us.

San Up S. A. guarantees that this product has been manufactured with high quality materials and fulfills all the standing safety standards, both national and international.

Please read the instructions very carefully before using it and use the machine as described here.

## Important Safeguards

---

1) Read the instruction manual carefully before using the unit and keep it for future consultation.

2) The "Teddy-Neb" is a nebulizer for the inhalation of medical aerosols to treat asthma, chronic bronchitis and other respiratory tract disorders. This nebulizer is suitable for all commonly used inhalation solutions. Use only the type, amount of medicine approved by your local health organization indicated by your doctor.

3) This nebulizer must be used only for the specific cause it was created for: the aerosoltherapy. Any other use will be considered unsuitable and thus dangerous; the manufacturer cannot be considered responsible for the damages caused by unsuitable, incorrect or irrational use, or if the unit has been used in electrical installations not in accordance with the standing safety standards.

4) The following precautions are especially important:

- The unit is not protected from liquid penetration.
- The nebulizer must not be touched with wet or soaked hands.
- Do not use in extreme weather conditions.
- Do not use the device in presence of inflammable anesthetic mix in contact with the air, oxygen or nitrogen protoxide.
- Do not block the air ventilation holes on left/right side of the compressor unit.
- Unplug the unit before cleaning or maintaining it.
- Disconnect the nebulizer when it is not in use.

Do not untwist the jet nebulizer when the unit is operating.

Do not leave medication inside the jet nebulizer when it is not in use.

Do not pour more than 10 ml in the jet nebulizer.

Use the unit on stable surfaces. Never use it while sleepy or drowsy.

Do not submerge the nebulizer in water.

Use the unit intermittently. Maximum operating period: 30 minutes. It must stand unused for 30 minutes.

-Children will use this unit only under adult supervision.


-The mask and the mouth piece are for personal use, and they should not be shared in order to avoid any infection or contagion.

-Stop the nebulization if you feel any indisposition.


-This unit have to be repaired only by authorized personal.

**NOTE: this unit employs a thermal limiter that will shut off the entire product should the motor overheat.**

## Features
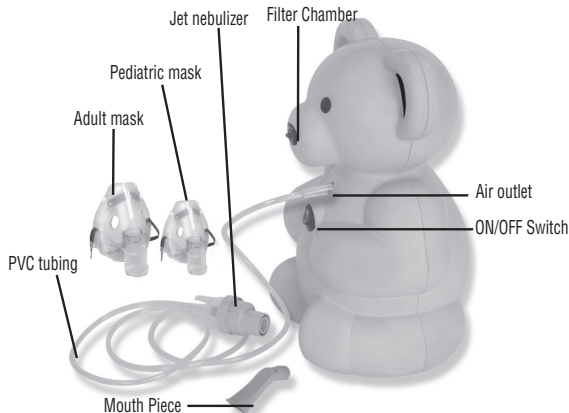
## Unit Assembly
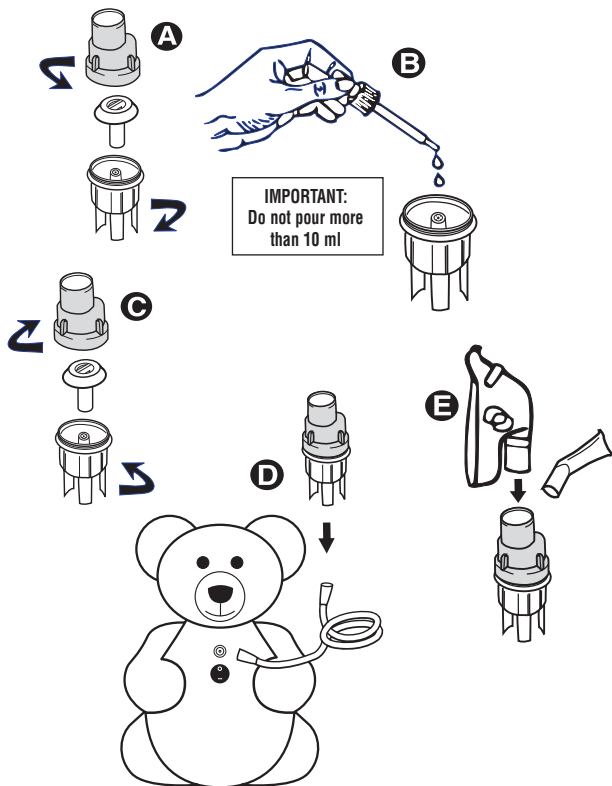
## Operating Instructions

---

- 1- Take the unit out of the box.
- 2- Remove the accessories.
- 3- Open the jet nebulizer (Fig. A) and pour the medication to nebulize into it. Do not pour more than 10 ml. (Fig. B), then twist both parts and close again. (Fig. C).
- 4- Connect the tubing to the jet nebulizer and the other end to the nebulizer outlet. (Fig. D).
- 5- Put the mask or the mouth piece to be used into the jet nebulizer (Fig. E).
- 6- Connect the unit to a suitable power receptacle.
- 7- Turn the unit by pressing the ON/OFF switch to start the treatment.
- 8- Turn the unit off and unplug it once the treatment has finished.
- 9- Clean the jet nebulizer and the devices as indicated in "Cleaning" section.
- 10- Store the jet nebulizer and the devices into the unit's lodging.

## Filter change

---

- 1- Filter should be changed every month or sooner if the filter turns completely gray in color.
- 2- Remove filter cap by grasping it firmly and pulling out the unit.
- 3- Change the dirty filter of the cap with small, pointed object.
- 4- Insert a clean filter. Extra filters are available.
- 5- Push the filter cap with new filter into position.

### Caution:

Do not use cotton or any other material. Do not wash or clean the filter. Use only San Up filters. Do not operate without a filter.

## Cleaning and Maintenance

---

### Important:

Before using the nebulizer for the first time or after each use the devices must be cleaned as follows: clean all the devices, except the tubing, with warm water and mild detergent solution. Rinse with warm tap water for 30 seconds to remove detergent. Put all the parts in a mixture of three parts water to one part white vinegar during 30 minutes into a container or a saucepan, or use a bacteria germicide disinfectant, which can be obtained in drugstores.

Clean your hands properly and take the parts out of the disinfectant solution, rinse them with warm tap water and let air dry on clean paper towels. Put them into a hermetic bag. Do not dry the nebulizer parts with a towel, this could contaminate them.

**NOTICE: It is recommended to clean all the parts after each use to avoid possible infection risk due to contaminated medication. Nebulize only products prescribed by the physician and follow his indications. Do not use balsamic or aromatic substances, they may cause the equipment and accessories to deteriorate.**

Use a dry cloth to clean the out part of the unit.

**The jet nebulizer and the devices must NOT be boiled.**

The jet nebulizer has to be replaced if it is stuck by medication or other substances.

Under no circumstances will the unit be opened. There is not any part inside it that could be repaired by non-qualified workmanship, the compressor does not require either lubrication or maintenance.

Do not use any volatile fluids, like benzene or alcohol, etc.

The life time of the accessories, using them according to the instructions of this manual is of 3 years.

## Spare parts provided

---

1 Jet nebulizer  
1 PVC tubing  
1 Mouth piece  
1 Air Filter (3)

1 Pediatric mask  
1 Adult mask  
1 Instruction Manual  
and Warranty

## Technical Specifications

Method of nebulization.....	Dry piston compressor
Operation voltage.....	<input type="checkbox"/> 110V~ <input type="checkbox"/> 120V~ <input type="checkbox"/> 220V~
Class.....	<input type="checkbox"/> 230V~ <input type="checkbox"/> 240V~ <input type="checkbox"/> 50Hz <input type="checkbox"/> 60Hz
	<input type="checkbox"/>
Power consumption.....	<input type="checkbox"/> 0.6A <input type="checkbox"/> 0.72A <input type="checkbox"/> 1.2A
Jet nebulizer maximum capacity.....	10 ml.
Sound noise level (at 1m. distance).....	55 dB
Mean Mass Aerodynamic Diameter (microns).....	less than 5 microns
Medication outlet.....	0.33 ml/min or greater*
Liter flow range.....	6-8,5 l/min
Approximate maximum pressure.....	225,56 kPa 32,73 PSI (2.3 kg/cm <sup>2</sup> )
Flow maximum .....	10 l/min
Weight.....	1.87 kg - 4,11lb
Unit dimensions.....	186 x 280 x 167 mm 7.32" x 11.02" x 6.57"
Operation Humidity and Temperature.....	5 to 40°C (41 to 104°F), 30% to 95%
Storage Humidity and Temperature.....	-10 to 40°C (19 to 104°F), 10% to 95%
Atmospheric pressure range for carriage and storage.....	50 - 106 kPa (500-1060 mbar)

\* value measured with 0.9% NaCl solution.

This specification may vary with the kind of the medication being used and it's viscosity.

	Pressure of work			Flow
	kg/cm <sup>2</sup>	PSI	kPa	l/min
Minimum	0,700	9,96	68,65	6
Maximum	0,950	13,52	93,16	8,5

• The specifications are subject to change without previous notice.

## Service Parts List

In case of requiring reposition, please ask for the parts mentioning the code number

Illustration	Description	Code number
	Jet nebulizer	1050
	PVC tubing	74346
	Pediatric mask	74207
	Adult mask	74208
	Mouth piece	74186
	Air filter	84005

## Symbol Key

On.....		Class II.....	
Off.....	○	Type B.....	
Alternative current.....	~	Read the instructions.....	

## Trouble Shooting Guide

If something unexpected happens during operation, please check and correct the following points:

Symptom	Possible cause	Action
When switching it on, it doesn't work.	1- Improperly plugged. 2- Electrical energy missing.	1- Check if the unit is properly connected 2- Verify if there is electrical energy.
The aerosol doesn't come out from the jet nebulizer.	1- Not enough medication in jet nebulizer. 2- Too much medication. 3- PVC tubing disconnected.	1- Add medication, maximum 10 ml. 2- Remove excess medication from cup. 3- Connect the PVC tubing.
Aerosol output reduced (low).	1- Air filter is clogged. 2- The jet nebulizer is dirty.	1- Replace for a clean air filter. 2- Replace for new one.

## Certifications

San Up S.A. meets the highest international quality standards and was awarded the following Certificates:

ISO 9001 Quality Systems (International)  
ISO 13485 Medical Equipment (European Community)

93/ 42/ CE Directive (European Community)

**A.N.M.A.T.**

(Argentina)

(Argentina)

San Up S.A. guarantees that all the supplies used in the nebulizer do not produce any risk due to substances released during its use.  
The nebulizer complies with the biocompatibility test according to ISO 10993.

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Tel.: (5411) 4754 4000 Fax: (5411) 4754 4004

A.N.M.A.T. P.M. 213-22

Technical Director Marta Camesella

**Made in Argentina**

CE Marking EU Authorized Representative:

Name: MEDES LIMITED

Address: 5 Beaumont Gate, Shenley Hill,

City: Radlett, Herts WD7 7AR,

Country: England

Tel.: +423-663-169205

Tel./Fax: +44 1923859810

## 1 Year Limited Warranty

---

San Up S.A. warrants to the original retail purchaser that this product is free from defects in material and workmanship for a period of 1 year from the date of purchase. This warranty does not cover failures resulting from misuse, abuse, alteration of the San Up Nebulizer System, use improper components or failure to comply with the operating instructions. This warranty does not include the Nebulizer System tubing or accessories.

THIS LIMITED WARRANTY IS IN LIEU OF ALL OTHER WARRANTIES, WHETHER ORAL, WRITTEN, EXPRESS OR IMPLIED. THERE ARE NO WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, EXCEPT FOR THOSE IMPLIED WARRANTIES ARISING BY OPERATION OF LAW WHICH ARE LIMITED IN DURATION TO THE 1 YEAR TERM OF THIS WARRANTY. SAN UP S.A. LIABILITY, WHETHER BASED ON CONTRACT, TORT, WARRANTY, STRICT LIABILITY, OR ANY OTHER THEORY, SHALL NOT EXCEED THE PRICE OF THE INDIVIDUAL SAN UP NEBULIZER SYSTEM WHICH IS THE BASIS OF THE CLAIM. IN NO EVENT SHALL SAN UP S.A. BE LIABLE FOR ANY LOSS OF PROFITS, FACILITIES, EQUIPMENT OR OTHER INCIDENTAL, CONSEQUENTIAL OR SPECIAL DAMAGES.

Some states do not allow limitations on how long an implied warranty lasts, or allow the exclusion or limitation of incidental, consequential, or special damages, so the above limitations or exclusions may not apply to you. This warranty gives you specific legal rights. You may also have other rights which vary from state to state.

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### Fabricado por San Up S.A.